

Bisturí y aspirador ultrasónicos Sonicmed



La tecnología ultrasónica de Sonicmed convierte la energía eléctrica en energía mecánica, haciendo vibrar las puntas a **velocidad ultrasónica a 40 kHz** de frecuencia. Esa potencia aumenta exponencialmente la eficiencia y precisión de operación minimizando el riesgo de dañar al paciente.

El sistema integrado de Sonicmed consta de **2 piezas de mano** separadas, ofreciendo un diseño compacto, cómodo y capaz de alcanzar la máxima potencia emitida por la consola.

- **Bisturí y fresas quirúrgicas**
- **Aspirador de tejido insano y dañado**



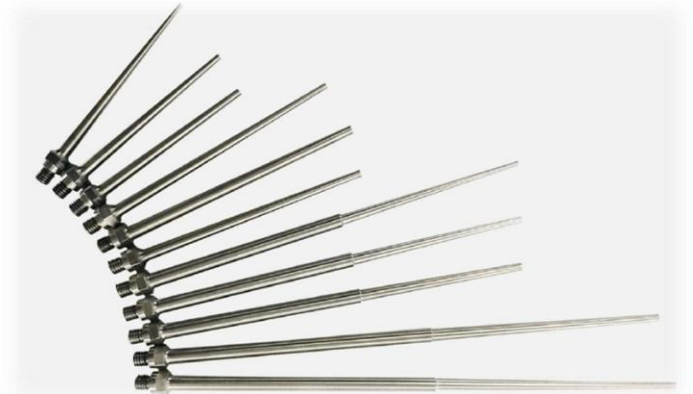
La tecnología selectiva del bisturí vibrando a velocidad ultrasónica está diseñada para cortar el hueso con facilidad sin dañar el tejido blando circundante, minimizando la pérdida de hueso y sangre del paciente.

El aspirador detecta, tritura y aspira el tejido dañado, agilizando cirugías tumorales sin pérdidas innecesarias de tejido sano.

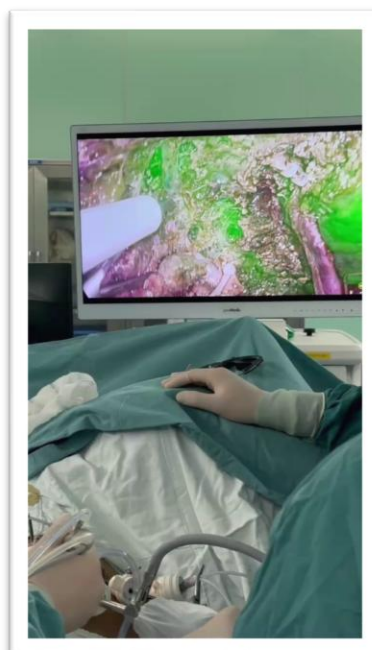
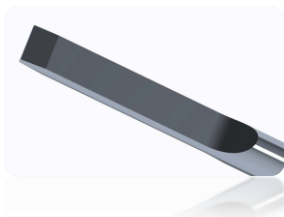


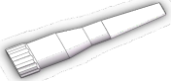
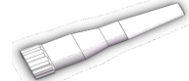
- **Amplia selección** de puntas de distintas formas, grosor y longitud.
- Vida útil media de **más de 40 minutos** de uso para las puntas de corte y **hasta 200 minutos** para las fresas.
- **Más de 40 puntas** distintas con la posibilidad de diseñar piezas personalizadas

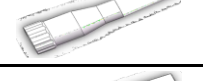
- Longitud entre **70 y más 300 mm** de longitud.
- Distintos diámetros de orificio de succión para una **aspiración más o menos agresiva**.
- **Vida útil** virtualmente ilimitada.

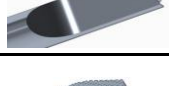
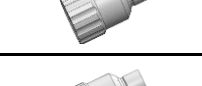
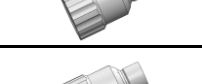
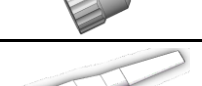


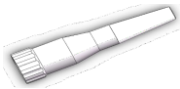
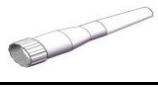
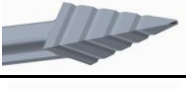
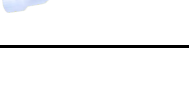
Además, las piezas de mayor longitud permiten el uso del sistema Sonicmed en cirugías asistidas con **endoscopia y mínimamente invasivas**.



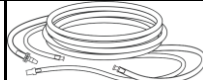
No.	Product type	Handpiece	Product code	Tip model	Tip Details	Product Picture	Description	Specification	Product code	Sleeve	Sleeve Picture
1	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402029	OT3CBN			Short	Length: 70mm Tip thickness: 3mm	407063	OSC40	
2	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402033	OT3CBL			Long	Length: 140mm Tip thickness: 3mm	407064	OSC8NL	
3	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402034	OT5CBL			Long	Length: 140mm, Tip thickness: 5mm	407064	OSC8NL	
4	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402026	OT3CNH			Short	Length: 70mm Tip thickness: 3.5mm	407066	OSD40	
5	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402028	OT2CNN			Short	Length: 70mm Tip thickness: 2.5mm	407063	OSC40	
6	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402032	OT2CNL			Long	Length: 140mm, Tip thickness: 2.5mm	407064	OSC8NL	
7	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402039	OT7CRT			Minimally invasive	Length: 300mm Tip thickness: 8mm	407065	OSD2NS	
8	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402027	OT5CBN			Short	Length: 70mm Tip thickness: 5mm	407063	OSC40	
9	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402035	OT3CNL			Long	Length: 140mm Tip thickness: 3.5mm	407067	OSD3NL	
10	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402038	OT2CNS			Minimally invasive	Length: 300mm Tip thickness: 3.65mm	407065	OSD2NS	
11	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402031	OT4CUN			Short	Length: 70mm Tip thickness: 4mm	407063	OSC40	

12	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402036	OT2CXS			Minimally invasive	Length: 310.9mm Tip thickness: 3.65mm	407065	OSD2NS	
13	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402197	OT3CXN12			Lengthen	Length: 192mm Tip thickness: 3mm		OS3CXN12	
14	Bone Drill Tip	40K Ultrasonic Bone Drill Handpiece	402021	OT4SPL			Long	Length: 140mm Tip thickness: 4mm	407067	OSD3NL	
15	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402001	OT8SNN			Short	Single-sided toothed, Length: 80mm Tip thickness: 0.8mm	407063	OSC40	
16	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402013	OT8SNL			Long	Length: 140mm Tip thickness: 0.8mm	407064	OSC8NL	
17	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402007	OT8SSN			Short	Length: 74.15mm Tip thickness: 0.8mm Tip Length: 15mm	407063	OSC40	
18	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402005	OT8SHN			Short	Length: 70mm Tip thickness: 0.8mm		OSC8HN	
19	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402002	OT8SFN			Short	Length: 70mm Tip thickness: 5mm		OSC8FN	
20	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402004	OT8SGN			Short	Length: 80mm Tip thickness: 0.8mm with a hole	407063	OSC40	
21	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402008	OT8SSH			Short	Length: 74.7mm Tip thickness: 0.8mm Tip Length: 20mm	407063	OSD40	
22	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402010	OT5SCN			Short	Length: 70mm Tip thickness: 0.8mm	407066	OSD40	
23	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402011	OT2SON			Short	Length: 70mm Tip thickness: 2.5mm	407063	OSC40	
24	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402012	OT8SNC			Short	Length: 70mm Tip thickness: 0.8mm	407065	OSD2NS	

25	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402014	OT8SFL			Long	Length: 140mm Tip thickness: 5mm	407064	OSC8NL	
26	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402015	OT8SHL			Long	Length: 140mm Tip thickness: 0.8mm	407064	OSC8NL	
27	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402016	OT3SWN			Short	Length: 70mm Tip thickness: 3mm	407063	OSC40	
28	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402017	OT3SWL			Long	Length: 140mm Tip thickness: 3mm	407064	OSC8NL	
29	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402018	OT8SSL			Long	Length: 140mm Tip thickness: 0.8mm	407064	OSC8NL	
30	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402019	OT2SOL			Long	Length: 140mm Tip thickness: 2.5mm	407064	OSC8NL	
31	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402022	OT8SSS			Minimally invasive	Length: 314mm Tip thickness: 3.65mm	407065	OSD2NS	
32	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402023	OT8SAS			Minimally invasive	Length: 251.7mm Tip thickness: 4mm		OSC8HN	
33	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402025	OT5SCS			Minimally invasive	Length: 320mm, Tip thickness: 1.2mm	407065	OSD2NS	
34	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece		OT8SIN11			Short	Teeth on both sides , Length: 74.15mm Tip thickness: 0.8mm	407063	OSC40	
35	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece		OT8SAN81			Short	Length: 73.1mm Tip thickness: 0.8mm with a hole	407063	OSC40	
36	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece		OT8SIN81			Short	Length: 71.15mm Tip thickness: 0.8mm with a hole	407063	OSC40	
37	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece-curved		OT8SIS11B			Short	Length: 74.15mm Tip thickness: 0.8mm	407063	OSC40	
38	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece-curved		OT8SIL12B			Long	Length: 138.5mm Tip thickness: 0.8mm	407064	OSC8NL	

39	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece		OT8SIN12			Short	Length: 74.7mm Tip thickness: 0.8mm	407063	OSC40	
40	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402003	OT4SFC			Short	Length: 70mm Tip thickness: 5mm		OSC8FN	
41	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402020	OT4STL			Long	Length: 140mm Tip thickness: 4mm	407064	OSC8NL	
42	Bone Cutting Tip	40K Ultrasonic Bone Cutting Handpiece	402009	OT5SSN			Short	Length: 70mm Tip thickness: 0.5mm	407066	OSD40	
43	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402044	AT34NM			Short	Length: 68.92mm Inner diameter: 2mm	407037	AA3IS3	
44	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402043	AT34NS			Short	Length: 71.4mm Inner diameter: 1.6mm	407041	AA34NS	
45	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402045	AT14NN			Short	Length: 75.8mm Inner diameter: 1mm	407037	AA3IS3	
46	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402048	AT34NN			Long	Length: 121.81mm Inner diameter: 2mm	407040	AA3IS4	
47	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402050	AT24NL			Long	Length: 125mm Inner diameter: 1.6mm	407040	AA3IS4	
48	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402049	AT14NL			Long	Length: 137mm Inner diameter: 1mm	407038	AA33NL	
49	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402193	AT34EN85			Long	Length: 123.95mm Inner diameter: 2mm		AS34EN85	
50	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402196	AT34NN55			Long	Length: 121.85mm Inner diameter: 2mm	407064	OSC8NL	

51	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402194	AT24NS74		Short	Length: 66.8mm Inner diameter: 1.66mm	407041	AA34NS	
52	Aspirator Tip	40K Ultrasonic Aspirator Handpiece	402195	AT54NN87		Short	Length: 64.5mm Inner diameter: 1.5mm	407041	AA34NS	
53	Aspirator Tip	40K Ultrasonic Aspirator Handpiece		AT34NT87		Minimally invasive	Length: 306.6mm Inner diameter: 2mm		AS34NT87	
54	Aspirator Tip	40K Ultrasonic Aspirator Handpiece		AT34NT88		Minimally invasive	Length: 367.8mm Inner diameter: 2mm		AS34NT88	
55	Aspirator Tip	40K Ultrasonic Aspirator Handpiece		AT34NT89		Minimally invasive	Length: 427.9mm Inner diameter: 2mm		AS34NT89	

56	Irrigation tube	40K Ultrasonic Bone Cutter & Aspirator Handpiece	SINGLE PIPELINE	
57	Irrigation and aspiration tube	40K Ultrasonic Bone Cutter & Aspirator Handpiece	DOUBLE PIPELINE	
58	Sterilization case	40K Ultrasonic Bone Cutter & Aspirator Handpiece	UDI06930934715524	

Distribuidor por:



BIOSpine. S.L.

Camino de los Cinamomos 283-5
33203 - Gijón (Asturias)
985195325

info@biospine.info